



Process Filtration From Pure to Sterile

(P)-PF-PES „W“

MAIN FEATURES & BENEFITS:

- Absolute ratings of 0,2 µm, 0,45 µm & 0,6 µm
- Outstanding flow rate
- Highly resistant materials
- Extremely low adsorption of proteins
- High thermal stability, permanently hydrophilic
- Approved for Food Contact Use
acc. to CFR Title 21 & 1935/2004/EC



INDUSTRIES:



- Bottled Water



- Soft Drinks



- Food & Beverage



- Biotech



- Chemical

Donaldson[®]
Ultrafilter

PRODUCT DESCRIPTION

The (P)-PF-PES „W“ filter element is an absolute rated, pleated high performance Polyethersulfone membrane filter. It provides the greatest assurance of filtration performance, stability and service life for sterile filtration and microbial stabilization.

The outstanding performance of the (P)-PF-PES „W“ filter element is based on its state-of-the-art filtration media. The Polyethersulfone membrane is inherently hydrophilic and distinguishes itself by having an asymmetrically designed pore structure. The pore size steadily decreases towards the centre of the medium resulting in a highly porous structure. This extremely durable design maintains consistent porosity and impurity retention throughout its service life without shedding or unloading contaminations.

All components meet the EU and USA requirements for Food Contact Use in accordance with **CFR (Code of Federal Regulations) Title 21 and 1935/2004/EC**. (P)-PF-PES „W“ has passed the USP Class VI tests for plastics. The filter element is manufactured in accordance with the manufacturing requirements, has no migration of filter media, is non-fibre releasing and is thermally welded without the use of binders or other chemical additives.

The absolute rated (P)-PF-PES „W“ membrane filter is designed and developed for the following applications:

Clarification and sterilization of all kinds of water

- Bottled Water
- Mineral Water
- Spring Water
- Table Water
- Potable Water

High quality filtration for a variety of ultrapure water requirements

- Deionized Water
- Chemically treated Water
- High temperature Water
- Process Water
- Ingredient Water

Sterile filtration of beverages

- Soft Drinks
- Beer
- Wine

PRODUCT SPECIFICATIONS

Product Specifications

Absolute Retention Rates

- 0,2 µm, 0,45 µm, 0,6 µm

Filtration Surface

- 0,72 m² per 250 mm element (10")

Maximum Differential Pressure

Operating temperature [°C / °F]	Differential pressure [bar / psi]
38 / 100	5,5 / 80
66 / 150	4,1 / 60
82 / 180	2,1 / 30

Cumulative Steaming Time*

- 121°C – 125°C (30 minutes) Saturated Steam (Forward Flow) up to 100 cycles

*Figures are based on lab tests to evaluate steaming resistance. Filter elements need to be checked in actual use. Contact Donaldson for recommended Autoclaving/Steaming procedures.

MATERIAL COMPLIANCE EU

The Donaldson (P)-PF-PES „W“ filter element meets the guideline for Food Contact Use as given in **European Regulation (EC) Number 1935/2004**. All polymeric components (Polypropylene, Polyethersulfone) meet the requirements of EU Directive 2002/72/EC relating to plastic materials and articles intended to come into contact with foodstuffs (excluding O-rings).

Migration tests have been carried out in simulant after flushing or in flow conditions.

For specific details on the O-rings, please contact your Donaldson Sales Engineer.

MATERIAL COMPLIANCE USA

All components of the (P)-PF-PES "W" filter element are FDA listed for food contact use in the **Code of Federal Regulations (CFR), Title 21**

Filter Materials	CFR Title
Membrane:	Polyethersulfone
Upstream Support:	Polypropylene
Downstream Support:	Polypropylene
Outer Guard:	Polypropylene
Core:	Polypropylene
End Caps:	Polypropylene
O-Rings:	EPDM
Alternatively:	Silicone
	Buna N
	PTFE over silicone
	PTFE over viton
Sealing Method:	Thermal Bonding

All products have been inspected and released by Quality Assurance as having met the following requirements:

- All filters are fabricated without the use of binders, adhesives, additives or surface-active agents.
- All filter components based on plastics are non-toxic and are certified bio-safe in accordance with current USP Class VI Tests for Plastic.
- All sterile filters are integrity tested to verify compliance with established quality and design specifications and to assure consistent and reliable performance.
- Bacterial endotoxin levels in aqueous extracts of (P)-PF-PES „W“ filter cartridges are less than 0,5 EU/ml, as determined using the limulus amebocyte lysate (LAL) test.


RETENTION RATES (ACCORDING TO HIMA CHALLENGE PER ASTM)

Filter Grade	Microorganism	LRV / cm ²
(P)-PF-PES „W“, 0,6 µm	<i>Saccharomyces cerevisiae</i>	> 7
	<i>Oenococcus oeni</i>	> 7
(P)-PF-PES „W“, 0,45 µm	<i>Saccharomyces cerevisiae</i>	> 7
	<i>Oenococcus oeni</i>	> 7
	<i>Serratia Marcescens</i>	> 7
(P)-PF-PES „W“, 0,2 µm	<i>Saccharomyces cerevisiae</i>	> 7
	<i>Oenococcus oeni</i>	> 7
	<i>Serratia Marcescens</i>	> 7
	<i>Pseudomonas diminuta</i>	> 7

INTEGRITY TESTING


Bubble Point Test		Diffusion Test / Forward Flow Test	
Filter Grade	Minimum Bubble Point [bar / psi]	Filter Grade	Maximum Diffusion Values [ml / min]
0,6 µm	1,24 bar / 18 psi	0,6 µm	15 ml/min @ 0,7 bar (10 psi)
0,45 µm	2,21 bar / 32 psi	0,45 µm	25 ml/min @ 1,7 bar (25 psi)
0,2 µm	3,03 bar / 44 psi	0,2 µm	30 ml/min @ 2,4 bar (35 psi)

FLOW CHARACTERISTICS

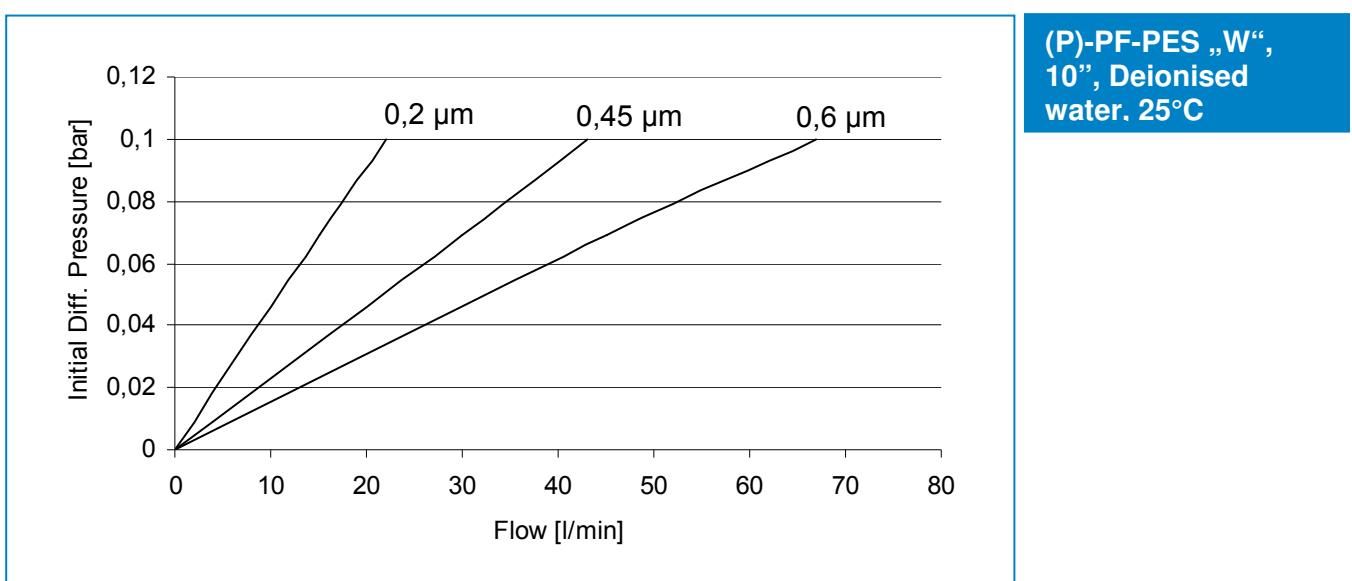

AVAILABLE END CAP CONFIGURATIONS**Dimensions (CODE 7 connection):**


Size	a		b		c	
	mm	inch	mm	inch	mm	Inch
10"	250	9,84	315	12,40	56,5	2,22
20"	500	19,68	565	22,24	56,5	2,22
30"	750	29,53	815	32,08	56,5	2,22

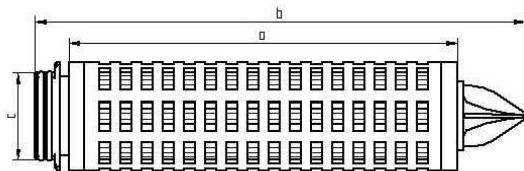

CODE 7: 2 x 226 o-rings, bayonet 2 locking tabs, locating fin.

Dimensions (uf plug connection):


Size	a		b		c	
	mm	inch	mm	inch	mm	Inch
10/30	254	10	270	10,63	61	2,40
20/30	510	20	526	20,63	61	2,40
30/30	764	30	780	30,63	61	2,40

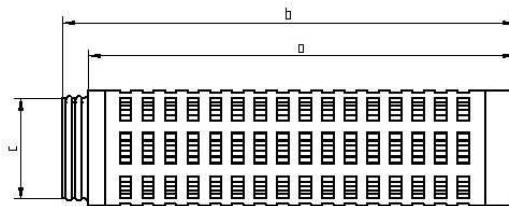

uf: plug connection with integrated reinforcement metal ring and flat end cap.

Dimensions (P9 connection):

Size	a		b		c	
	mm	inch	mm	inch	mm	Inch
10"	250	9,84	320	12,59	44	1,73
20"	500	19,68	570	22,44	44	1,73
30"	750	29,53	820	32,28	44	1,73

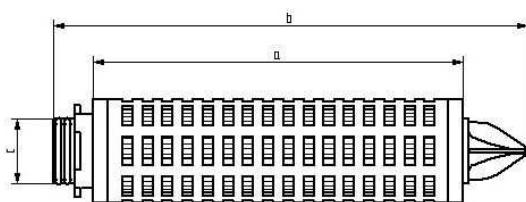

P9: 2 x 222 o-rings, bayonet 3 locking tabs, locating fin.

Other end cap configurations on request.

Technical alterations reserved 04/2009

- **Integrity test to be done by Bubble Point or Forward Flow Test.**
- **For information on test equipment or test services, please contact your Donaldson Sales Engineer**