



Donaldson  
FILTRATION SOLUTIONS

# Process Filtration From Pure to Sterile

## PF-PES "B"

### MAIN FEATURES & BENEFITS:

- Absolute ratings of 0,2  $\mu\text{m}$  & 0,45  $\mu\text{m}$
- Outstanding flow rate
- Highly resistant materials
- Extremely low absorption of proteins
- High thermal and hydrolytic stability
- Approved for Food Contact Use acc. to CFR Title 21 & 1935/2004/EC



### INDUSTRIES:



- Breweries



- Wineries



- Mineral Water



- Soft Drinks



- Chemical



- Pharmaceutical

Donaldson®  
**Ultrafilter**

## PRODUCT DESCRIPTION

---

The PF-PES „B“ filter element is an absolute rated, pleated high performance Polyethersulfone membrane filter. It provides the greatest assurance of filtration performance, stability and service life for sterile filtration and microbial stabilization.

---

The outstanding performance of the PF-PES „B“ filter element is based on its state-of-the-art filtration media. The Polyethersulfone membrane is inherently hydrophilic and distinguishes itself by having an asymmetrically designed pore structure. The pore size steadily decreases towards the centre of the medium resulting in a highly porous structure. This extremely durable design maintains consistent porosity and impurity retention throughout its service life without shedding or unloading contaminations.

---

This extremely durable design maintains consistent porosity and impurity retention throughout its service life without shedding or unloading contaminations.

---

All components meet the EU and USA requirements for Food Contact Use in accordance with **CFR (Code of Federal Regulations) Title 21** and **1935/2004/EC**.

PF-PES „B“ has passed the USP Class VI tests for plastics. The filter element is manufactured in accordance with the manufacturing requirements, has no migration of filter media, is non-fibre releasing and is thermally welded without the use of binders or other chemical additives.

**The absolute rated PF-PES „B“ membrane filter is designed and developed for the following applications:**

### Clarification and cold sterilization of beverages like:

- Beer
- Beer Mix
- Wine
- Wine Coolers

### Clarification and final filtration of

- Deionized Water
- Chemically treated Water
- High temperature Water
- Process Water
- Ingredient Water
- Soft Drinks
- Bottled Water

## PRODUCT SPECIFICATIONS

## Product Specifications

## Absolute Retention Rates

- 0,2 µm, 0,45 µm

## Filtration Surface

- 0,72 m<sup>2</sup> per 250 mm element (10")

## Maximum Differential Pressure

Operating temperature [°C / °F]	Differential pressure [bar / psi]
38 / 100	5,5 / 80
66 / 150	4,1 / 60
82 / 180	2,1 / 30

## Cumulative Steaming Time\*

- 121°C – 125°C (30 minutes) Saturated Steam (Forward Flow) up to 100 cycles\*

*\*Figures are based on lab tests to evaluate steaming resistance. Filter elements need to be checked in actual use. Contact Donaldson for recommended Autoclaving/Steaming*

## MATERIAL COMPLIANCE EU

The Donaldson PF-PES „B“ filter element meets the guideline for Food Contact Use as given in **European Regulation (EC) Number 1935/2004**. All polymeric components (Polypropylene, Polyethersulfone) meet the requirements of EU Directive 2002/72/EC relating to plastic materials and articles intended to come into contact with foodstuffs (excluding O-rings).

Migration tests have been carried out in simulants after flushing or in flow conditions.

For specific details on the O-rings, please contact your Donaldson Sales Engineer.

**MATERIAL COMPLIANCE USA**

All components of the PF-PES „B“ filter element are FDA listed for food contact use in the **Code of Federal Regulations (CFR), Title 21**

Filter Materials	CFR Title
<b>Membrane:</b>	Polyethersulfone 177.2240
<b>Upstream Support:</b>	Polypropylene 177.1520
<b>Downstream Support:</b>	Polypropylene 177.1520
<b>Outer Guard:</b>	Polypropylene 177.1520
<b>Core:</b>	Polypropylene 177.1520
<b>End Caps:</b>	Polypropylene 177.1520
<b>O-Rings:</b>	EPDM 177.2600
<b>Alternatively:</b>	Silicone 177.2600
	Buna N 177.2600
	PTFE over silicone 177.1550
	PTFE over viton 177.1550
<b>Sealing Method:</b>	Thermal Bonding

**All products have been inspected and released by Quality Assurance as having met the following requirements:**

- All filters are fabricated without the use of binders, adhesives, additives or surface-active agents.
- Bacterial endotoxin levels in aqueous extracts of PF-PES „B“ filter elements are less than 0,5 EU/ml, as determined using the limulus amebocyte lysate (LAL) test.
- All filter components based on plastics are non-toxic and are certified bio-safe in accordance with current USP Class VI Tests for Plastic.
- All sterile filters are integrity tested to verify compliance with established quality and design specifications and to assure consistent and reliable performance.

RETENTION RATES (ACCORDING TO HIMA CHALLENGE PER ASTM)

Filter Grade	Microorganism	LRV / cm <sup>2</sup>
PF-PES „B“, 0,45 µm	Saccharomyces cerevisiae	> 7
	Oenococcus oeni	> 7
PF-PES „B“, 0,2 µm	Saccharomyces cerevisiae	> 7
	Oenococcus oeni	> 7
	Serratia Marcescens	> 7

INTEGRITY TESTING

Bubble Point Test


Filter Grade	Minimum Bubble Point [bar / psi]
0,45 µm	1,38 bar / 20 psi
0,2 µm	2,3 bar / 44 psi

Diffusion Test / Forward Flow Test

Filter Grade	Maximum Diffusion Values [ml / min]
0,45 µm	15 ml / min @ 0,7 bar (10 psi)
0,2 µm	25 ml / min @ 1,7 bar (35 psi)

FLOW CHARACTERISTICS


PF-PES „B“, 10”, Deionized water, 25°C

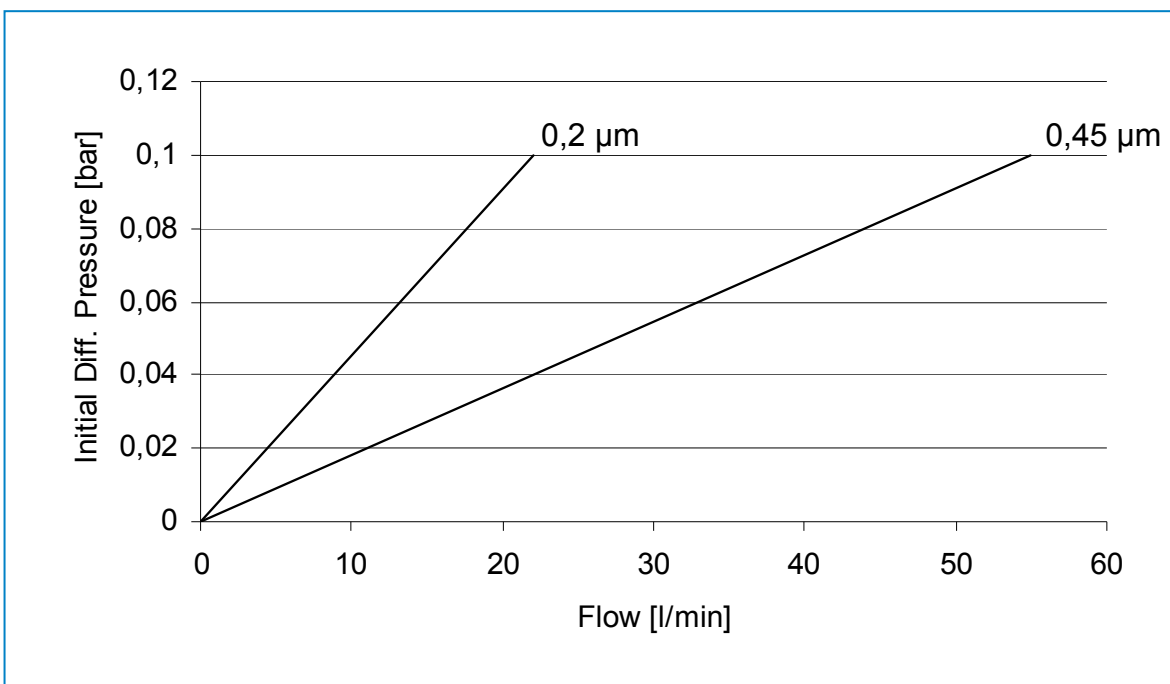

AVAILABLE END CAP CONFIGURATIONS

Dimensions (Code 7 connection)


CODE 7						
Size	a		b		c	
	mm	inch	mm	inch	mm	inch
10"	250	9,84	315	12,40	56,5	2,22
20"	500	19,68	565	22,24	56,5	2,22
30"	750	29,53	815	32,08	56,5	2,22

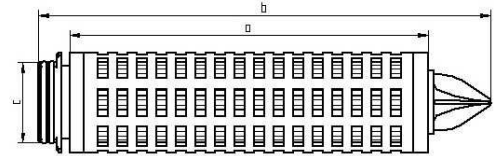

CODE 7: 2 x 226 o-rings, bayonet 2 locking tabs, locating fin

Dimensions (P9 connection)

Code 9						
Size	a		b		c	
	mm	inch	mm	inch	mm	inch
10"	250	9,84	320	12,59	44	1,73
20"	500	19,68	570	22,44	44	1,73
30"	750	29,53	820	32,28	44	1,73

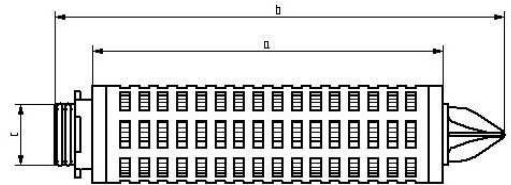

P9: 2 x 222 o-rings, bayonet 3 locking tabs, locating fin.

Other end cap configurations on request.

Technical alterations reserved 04/2009

- Integrity test of this element to be done by Bubble Point or Forward Flow Test.
- For information on test equipment or test services, please contact your Donaldson Sales Engineer

(Rev02 – 07/10)

